

# Apply the Force to your Consultations!

Dental Clinicians can use the knowledge of the force to identify potential restoration limitations.

Christopher Sitaras, M. Eng.;  
Jocelyn Vignola, D.D



**Innobyte™**  
by Kube Innovation



KUBE INNOVATION



# Background

Bite force is a key indicator of a healthy masticatory system, in that a deficiency or surplus is implicated in multiple disorders and dental complications. Amongst others, these disorders and complications are Temporomandibular Disorder (TMD),<sup>1</sup> bruxism,<sup>2</sup> restoration failure,<sup>3</sup> bone resorption,<sup>4</sup> and neurological diseases.<sup>5</sup> Of particular interest to clinicians is oral pain and the failure of dental restorations, the latter of which is chiefly caused by occlusal loading conditions, including bite force.<sup>6</sup>

As many as 35% of restorations fail,<sup>6-8</sup> which is a clinician's worst nightmare<sup>9</sup> because most of the time, the clinicians redo the work pro bono. The estimated costs of reparation during its lifetime can be as much as 12,000\$ per restoration that cannot always be recouped.<sup>4</sup> Additionally, there is an innumerable cost to a clinician's reputation, as the majority of patients suffering from failures consult a different practitioner.<sup>10</sup>

Currently, clinicians rely on a qualitative evaluation of a patient's bite force in order to provide a treatment plan. However, this measure is subjective and does not adequately characterize the magnitude of force outputted by the patient. A clear understanding of this value will assist the dental clinician in designing a personalized solution that adequately considers the imposed limitations of these occlusal loads.



# Understanding Bite Force

The forces from the masticatory muscles transfer their strength through the dentition during occlusion, and while some vertical load is applied, the principal action is axially.<sup>11,12</sup>

This load acts on the occlusal contact area and has been shown to be up to three times greater in the posterior jaw versus the anterior. When occlusal contact area is decreased (i.e.: missing teeth or malocclusion) there can be up to 20% decrease in bite force.<sup>13,14</sup> Bite force for all ages and genders vary greatly between 10 to 2000 N,<sup>15</sup> however, men are consistently found to have up to twice the levels found in women.<sup>16,17</sup> The average total bite force has been found to be 800 N and 600 N for adult male and female patients, respectively, decreasing gradually through to old age.<sup>13,18,23</sup> Also, single tooth average bite forces have been estimated to being 300 N for males and 250 N for females,<sup>12,19</sup> with maxima of as much as 1600 N.<sup>20</sup>

The age of dentate patients was closely tied to the maximal bite force, with a lower force recorded in children and the elderly, and a maximal force in early adulthood tapering off with increasing age, attributed to the decline in musculature with age.<sup>16,21</sup> It has been shown that the molar regions exhibit higher forces regardless of gender and age.<sup>16,22</sup>

These regions experience the highest force values, with as much as 67% of the total bite force in the posterior of the oral cavity,<sup>12,14</sup> and up to 58% of the whole arch bite force on the preferred chewing side.<sup>23</sup> Despite loss of dentition, the region could exert similar forces upon restorations due to comparable muscle function to dentate patients.<sup>24</sup> Despite varied levels of bite force in the general population, a clear understanding of the service conditions must be considered during the placement of the dental restoration.<sup>25</sup>

Crucial among these service conditions is the maximum bite force, which dictates the capacity a dental restoration should be able to withstand during regular use, to ultimately reduce the risk of failures.



# Bite Force and Failure in Dental Restorations

For clinicians to provide dental restorative treatments to patients, the bite force must be known in order to avoid restoration failures, which can be costly to repair and negatively impact reputations.<sup>4</sup> Amongst the etiological causes of restoration failures the highest remains mechanical fracture.<sup>3,26-30</sup> The cause of any material fracture is excessive stress on the component due to unfavorable loading conditions.

Under in vitro static loading conditions, a standard crown restoration fractures at a load of approximately 900 to 2400 N depending on the material used, type of preparation and many other factors.<sup>31,32</sup> Implants fracture at a load of approximately 150 to well over 1000 N.<sup>33</sup>

Bite force values of this magnitude are uncommon on a single tooth, however, under in vitro cyclical forces, the fatigue corrected loading conditions required to cause failure could be as low as 400 N,<sup>34,35</sup> which is within the range of most patients.<sup>12</sup> The cyclical nature of chewing implies that an occlusal load is applied to restorations as many as 3500 times daily,<sup>36</sup> which is an important factor in fatigue failure of material due to crack propagation.<sup>28,31,37,38</sup>

This is exacerbated by the environment in vivo which is both moist and cyclically heated, two parameters that can decrease the fatigue strength of dental materials by as much as 25%.<sup>39</sup> [The translation to years of service for some dental restoration, as a result of this environment, can be less than 1 year.](#)<sup>40</sup>

Thus, when force is being applied to the restoration that exceeds the materials' fatigue strengths, failure can occur earlier than its prescribed service life.

This means bite force must be considered when choosing the correct treatment plan for the patient in order to prevent dental restoration failures by assuring the choice of material can withstand the service conditions during its prescribed lifetime.

*Though the frequency of reported failures among dental restorations remains high,<sup>6</sup> compounding effects can increase the total failure rate to as much as 34%,<sup>9,41</sup>*

A survey of dentists in the Dental Practice-Base Research Network<sup>10</sup> have shown that [more than 50% of all failures occurred in the posterior region, where the largest of forces are applied, and dental amalgams fail in more than 50% of restorative failures, which is considered amongst the weaker dental restorative materials.](#)<sup>42</sup>

These results have been demonstrated separately using a similarly large cohort, which objectively shows that failures are region specific and material specific.<sup>43</sup> This strongly suggests that material and loading conditions play a crucial role in the service life of a restorative treatment plan.

These failure rates are based on curated reviews of clinical cases, while in practice, where most variables cannot be controlled, failure can occur more often.<sup>44</sup> While failure can be multifactorial, an implant<sup>35</sup> or crown<sup>45</sup> material whose fatigue strength is exceeded from an excessive bite force will not have a high survival rate.

Although an understanding of the restoration geometry is important for its application, a lack of understanding in material strength and forces increases the risk of a restoration's failure. To draw a professional parallel, engineers could not provide safe structural solutions for a bridge without accounting for the loading conditions and the materials assuming the load.

Currently, a dental clinician must not only choose the correct materials from a plethora of choices,<sup>46</sup> they must also assume the patient's loading conditions using uniquely qualitative clinical judgment, with very little quantitative information.

Though qualitative evaluation can be crucial in the installation of successful dental restorations, the high instances of failure may be avoided with a quantitative understanding of the patient's masticatory load.<sup>3,47</sup>



A failed dental restoration costs time and money, with as much as 12,000\$ in lifetime repair costs,<sup>4</sup> and approximately 30 minutes per repair of potentially unbillable work.<sup>48</sup> Most importantly, the reputation of the clinician is negatively impacted, which usually cannot be recouped, as the majority of reparation of dental restorations are effected by different clinicians.<sup>10</sup>

## A Case for Night Guards

The average patient applies approximately 800 N voluntarily on a single bite,<sup>12,23</sup> with potentially larger forces applied involuntarily during sleep parafunction events (eg: bruxism).<sup>2</sup> As some patients may exceed these average voluntary bite force values, so too does their bite force increase during bruxism events.<sup>49</sup>

As the magnitude and frequency of force is increased, there is an increased risk to dental restorations. In a review of the survival rate for dental implants in people with bruxism, it was found that there was two to four time odd ratio of implant failure in bruxers than those without bruxism.<sup>50</sup>

As the cost of dental restorations in time and money is significant, for patient and clinician,<sup>51,52</sup> effective protection is a welcomed addition to any dental restoration treatment plan, especially in instances of elevated bite force value.



## Solution: Personalized Bite Force Assessment



*The Innobyte™ is the first total bite force measurement device on the market. With the Innobyte's patent pending technology, the clinician is able to measure the whole arch loading from 0 to 3000 Newtons.*

Using the pressure of an incompressible fluid, the Innobyte is able to convert this value to an easy-to-read, easy-to-interpret, digital read-out of the calibrated force in Newtons (1 Newton is equal to 0.225-pound force). Due to the soft nature of the medical grade silicone of the bite surface, patients do not risk pain or damage to teeth or gingiva, and are able to apply their maximal load.<sup>53</sup>

The intermolar separation of the Innobyte's Mouthpiece is specially design in order to minimize elevator muscle bias, where less muscle fibers are recruited during large jaw separations,<sup>54,55</sup> as well as simulating the clinically relevant average food morsel size. The Innobyte's Mouthpiece is a one-size-fits-all flexible unit, so there is no need to worry about the fit with each patient.

The Innobyte is simplistic in design and operation, with no complex accessories and components, with an intuitive user interface. It has never been so easy to assess the maximum voluntary bite force of a patient.



# Accuracy, Precision and Reproducibility

The Innobyte associates the pressure from proprietary calibration system in order to provide precise and reproducible results. This is achieved by using calibrated system sensors and load cells that are based on ISO 17025 precision standards with less than 1% variation, for a compounded error of less than 5% (see Figure 1). This translates to a maximum variation of 20 Newtons or 5-pound force, or less than the force required to tap a key on a keyboard, which is quite impressive when the average human bites is well over 500 Newtons of force. Additionally, the Mouthpiece has four carefully designed guides for precise Mouthpiece placement that allows the clinician to obtain reproducible results, every time (see Figure 2). Differences in the measured bite force can be confidently associated with a successful treatment plan.

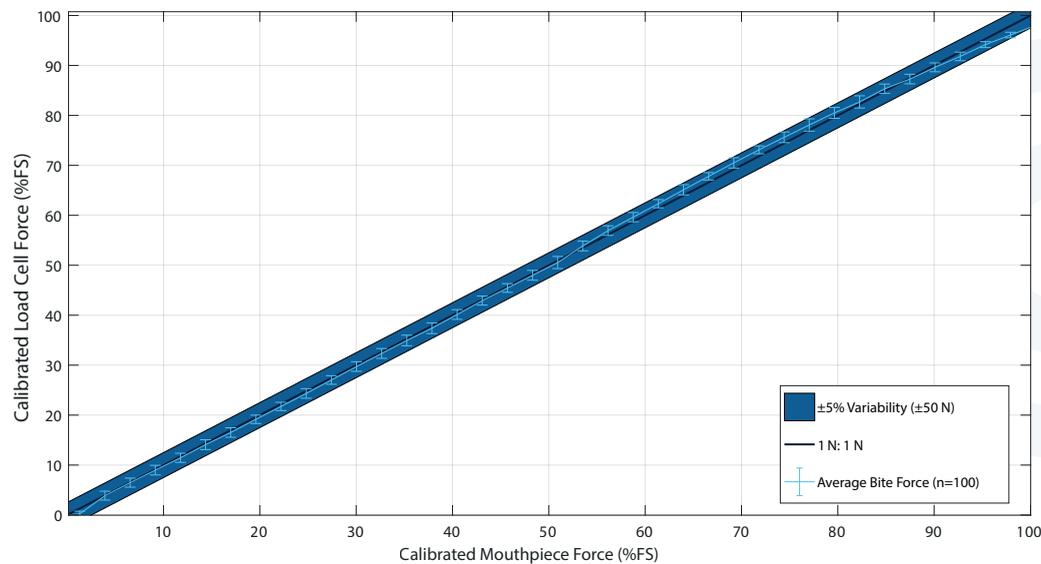


Figure 1. Calibration Curves for Innobyte™ Mouthpieces from 0-100% of Full Scale (2000N).

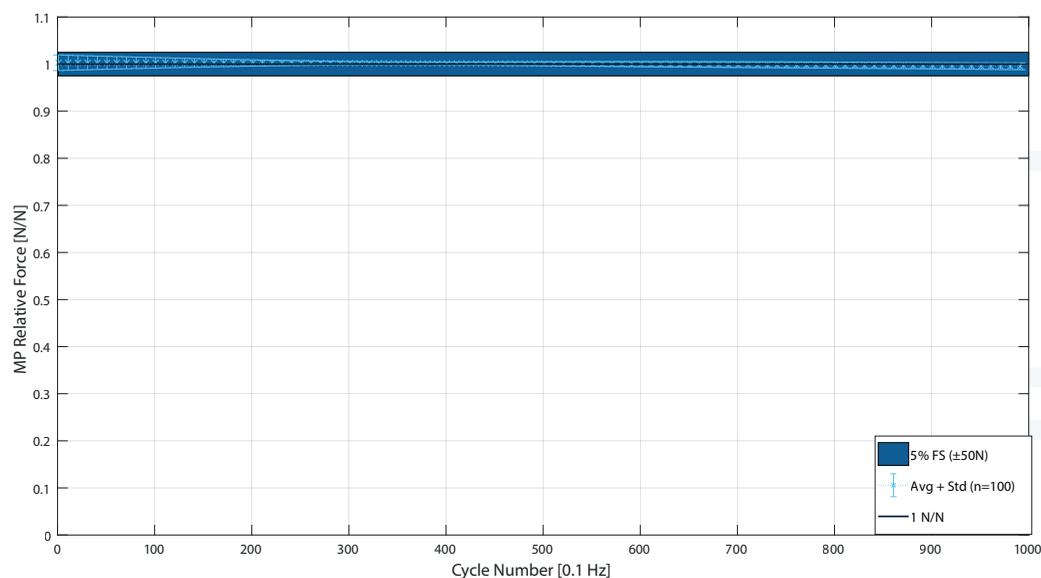


Figure 2. Repeatability of Bite Forces over entire lifecycle of Innobyte™ Mouthpiece (1000 cycles).



# Conclusion

The inclusion of bite force in clinical practice is an invaluable tool, which can help convince the patient that the proposed treatment plan is the correct course. Prior to the Innobyte, clinicians did not have all the tools to make a data driven decision for treatment plans. And due to its operating principles and ease of use, the Innobyte™ adds valuable objective information related to masticatory function in everyday clinical practice. Furthermore, bite force measurement can be easily integrated in everyday practice with the Innobyte™ to help in the detection and monitoring of the oral health of the clinician's patients.



## HIGHLY PRECISE

This high precision measuring device is over 95% accurate



## EDUCATIONAL

Use data to educate your patients



## USER FRIENDLY

Only press one button to get a precise bite force measure



# References

1. Clark GT, Beemsterboer PL, Jacobson R. The effect of sustained submaximal clenching on maximum bite force in myofascial pain dysfunction patients. *J Oral Rehabil.* 1984;11(4):387-391. doi:10.1111/j.1365-2842.1984.tb00590.x
2. Nishigawa K, Bando E, Nakano M. Quantitative study of bite force during sleep associated bruxism. *J Oral Rehabil.* 2001;28(5):485-491. doi:10.1046/j.1365-2842.2001.00692.x
3. Flanagan D. Bite force and dental implant treatment: a short review. *Med Devices.* 2017;10:141-148. doi:10.2147/M-DER.S130314
4. Misch CE. Age-Related Tooth Loss. In: *Text Book on Dental Implant Prosthetics.* 3rd ed. Elsevier; 2014:1-17.
5. Teixeira FB, Pereira Fernandes L de M, Tavares Noronha PA, et al. Masticatory deficiency as a risk factor for cognitive dysfunction. *Int J Med Sci.* 2014;11(2):209-214. doi:10.7150/ijms.6801
6. Oh SL, Shiau HJ, Reynolds MA. Survival of dental implants at sites after implant failure: A systematic review. *J Prosthet Dent.* 2020;123(1):54-60. doi:10.1016/j.prosdent.2018.11.007
7. Moraschini V, Barboza E d. SP. Success of dental implants in smokers and non-smokers: a systematic review and meta-analysis. *Int J Oral Maxillofac Surg.* 2016;45(2):205-215. doi:10.1016/j.ijom.2015.08.996
8. Chen H, Liu N, Xu X, Qu X, Lu E. Smoking, Radiotherapy, Diabetes and Osteoporosis as Risk Factors for Dental Implant Failure: A Meta-Analysis. *PLoS One.* 2013;8(8):e71955. doi:10.1371/journal.pone.0071955
9. Kate MA, Palaskar S, Kapoor P. Implant failure: A dentist's nightmare. *J Dent Implant.* 2016;6(2):51. doi:10.4103/0974-6781.202154
10. GROUP DC. DPBRN 5 : Reasons for replacement or repair of dental restorations Results: Overall and by Region. *J Am Dent Assoc.* Published online 2011. doi:https://doi.org/10.14219/jada.archive.2012.0238
11. Lundgren D, Laurell L. Occlusal forces in prosthetically restored dentitions: a methodological study. *J Oral Rehabil.* 1984;11(1):29-37. doi:https://doi.org/10.1111/j.1365-2842.1984.tb00549.x
12. Hattori Y, Satoh C, Kunieda T, Endoh R, Hisamatsu H, Watanabe M. Bite forces and their resultants during forceful intercuspal clenching in humans. *J Biomech.* 2009;42(10):1533-1538. doi:10.1016/j.jbiomech.2009.03.040
13. Bakke M. Bite Force and Occlusion. *Semin Orthod.* 2006;12(2):120-126. doi:10.1053/j.sodo.2006.01.005
14. Bakke M. Mandibular elevator muscles: physiology, action, and effect of dental occlusion. *Eur J Oral Sci.* 1993;101(5):314-331. doi:10.1111/j.1600-0722.1993.tb01127.x
15. Pruims GJ, de Jongh HJ, Ten Bosch JJ. Forces Acting on the Mandible During Bilateral Static Bite At Different Bite Force Levels. *J Biomech.* 1978;13(9):755-763. doi:https://doi.org/10.1016/0021-9290(80)90237-7
16. Helkimo E, Carlsson GE, Helkimo M. Bite Force and State Dentition. *Acta Odont Scand.* 1976;35(6):297-303. doi:10.3109/00016357709064128
17. Zivko-Babic J, Panduric J, Jerolimov V, Mioc M, Pizeta I, Jakovac M. Bite force in subjects with complete dentition. *Coll Antropol.* 2002;26(1):293-302.
18. Van Der Bilt A, Tekamp A, Van Der Glas H, Abbink J. Bite force and electromyography during maximum unilateral and bilateral clenching. *Eur J Oral Sci.* 2008;116(3):217-222. doi:10.1111/j.1600-0722.2008.00531.x
19. Takaki P, Vieira M, Bommarito S. Maximum bite force analysis in different age groups. *Int Arch Otorhinolaryngol.* 2014;18(3):272-276. doi:10.1055/s-0034-1374647
20. Jansen van Vuuren L, Broadbent JM, Duncan WJ, Waddell JN. Maximum voluntary bite force, occlusal contact points and associated stresses on posterior teeth. *J R Soc New Zeal.* 2020;50(1):132-143. doi:10.1080/03036758.2019.1691612
21. Bakke M, Holm B, Jensen BL, Michler L, Möller E. Unilateral, isometric bite force in 8-68-year-old women and men related to occlusal factors. *Eur J Oral Sci.* 1990;98(2):149-158. doi:10.1111/j.1600-0722.1990.tb00954.x
22. Waltimo A, Könönen M. A novel bite force recorder and maximal isometric bite force values for healthy young adults. *Eur J Oral Sci.* 1993;101(3):171-175. doi:10.1111/j.1600-0722.1993.tb01658.x
23. Khan SIR, Rao D, Ramachandran A, Ashok BV. Comparison of bite force on the dominant and nondominant sides of patients with habitual unilateral chewing: A pilot study. *Gen Dent.* 2020;68(2):60-63.
24. Haraldson T, Karlsson U, Carlsson GE. Bite force and oral function in complete denture wearers. *J Oral Rehabil.* 1979;6(1):41-48. doi:10.1111/j.1365-2842.1979.tb00403.x
25. Flanagan D. Diet and implant complications. *J Oral Implantol.* 2016;42(3):305-310. doi:10.1563/aid-joi-D-15-00127
26. Goiato MC, Haddad MF, Filho HG, et al. Dental Implant Fractures – Aetiology , Treatment and Case Report. *J Clin Diagnostic Res.* 2014;8(3):300-304. doi:10.7860/JCDR/2014/8074.4158
27. Stoichkov B, Kirov D. Analysis of the causes of dental implant fracture: A retrospective clinical study. *Quintessence Int (Berl).* 2018;49(4):279-286. doi:10.3290/j.qi.a39846
28. Zhang Y, Sailer I, Lawn BR. Fatigue of dental ceramics. *J Dent.* 2013;41(12):1135-1147. doi:doi:10.1016/j.jdent.2013.10.007.
29. Sakka S, Baroudi K, Nassani MZ. Factors associated with early and late failure of dental implants. *J Investig Clin Dent.* 2012;3(4):258-261. doi:10.1111/j.2041-1626.2012.00162.x
30. Lipski T, Kijak E, Witaszek K, Gańska D, Margielewicz J. Application of electronic dental dynamometer in biomechanics. *J Meas Eng.* 2015;3(1):17-22.



# References

31. Choi S, Yoon HI, Park EJ. Load-bearing capacity of various CAD/CAM monolithic molar crowns under recommended occlusal thickness and reduced occlusal thickness conditions. *J Adv Prosthodont.* 2017;9(6):423-431. doi:10.4047/jap.2017.9.6.423
32. Preis V, Hahnel S, Behr M, Bein L, Rosentritt M. In-vitro fatigue and fracture testing of CAD/CAM-materials in implant-supported molar crowns. *Dent Mater.* 2017;33(4):427-433. doi:10.1016/j.dental.2017.01.003
33. Bourauel C, Mundt T, Stark H, Heinemann F. Biomechanics and load resistance of small-diameter and mini dental implants: a review of literature. *Biomed Eng (NY).* 2013;59(1):1-5. doi:10.1515/bmt-2013-0092
34. Schultheis S, Strub JR, Gerds TA, Guess PC. Monolithic and bi-layer CAD/CAM lithium-disilicate versus metal-ceramic fixed dental prostheses: Comparison of fracture loads and failure modes after fatigue. *Clin Oral Investig.* 2013;17(5):1407-1413. doi:10.1007/s00784-012-0830-1
35. Kumararama SS, Chowdhary R. Selection of Dental Implants Based on Masticatory Load of the Patient: A Novel Approach. *Indian J Dent Res.* 2017;28(3):309-313. doi:10.4103/ijdr.IJDR
36. Fukuda H, Saito T, Mizuta M, et al. Chewing number is related to incremental increases in body weight from 20 years of age in Japanese middle-aged adults. *Gerodontology.* 2013;30(3):214-219. doi:10.1111/j.1741-2358.2012.00666.x
37. Zhang Y, Mai Z, Barani A, Bush M, Lawn BR. Fracture-resistant monolithic dental crowns. *Dent Mater.* 2016;32(3):442-449. doi:10.1016/j.dental.2015.12.010.FRACTURE-RESISTANT
38. Shahmoradi M, Wan B, Zhang Z, Wilson T, Swain M, Li Q. Monolithic crowns fracture analysis: The effect of material properties, cusp angle and crown thickness. *Dent Mater.* 2020;36(8):1038-1051. doi:10.1016/j.dental.2020.04.022
39. Rosentritt M, Preis V, Behr M, Strasser T. Fatigue and wear behaviour of zirconia materials. *J Mech Behav Biomed Mater.* 2020;110:103970. doi:10.1016/j.jmbbm.2020.103970
40. Moy PK, Medina D, Shetty V, Aghaloo TL. Dental implant failure rates and associated risk factors. *Int J Oral Maxillofac Implants.* 2005;20(4):569-577. <http://www.ncbi.nlm.nih.gov/pubmed/16161741>
41. Chappuis V, Avila-Ortiz G, Araújo MG, Monje A. Medication-related dental implant failure: Systematic review and meta-analysis. *Clin Oral Implants Res.* 2018;29:55-68. doi:10.1111/clr.13137
42. Chun KJ, Lee JY. Comparative study of mechanical properties of dental restorative materials and dental hard tissues in compressive loads. *J Dent Biomech.* 2014;5:1-6. doi:10.1177/1758736014555246
43. Van Nieuwenhuysen JP, D'Hoore W, Carvalho J, Qvist V. Long-term evaluation of extensive restorations in permanent teeth. *J Dent.* 2003;31(6):395-405. doi:10.1016/S0300-5712(03)00084-8
44. Da Silva JD, Kazimiroff J, Papas A, et al. Outcomes of implants and restorations placed in general dental practices: A retrospective study by the Practitioners Engaged in Applied Research and Learning (PEARL) Network. *J Am Dent Assoc.* 2014;145(7):704-713. doi:10.14219/jada.2014.27
45. Rosentritt M, Schumann F, Krifka S, Preis V. Influence of zirconia and lithium disilicate tooth-or implant-supported crowns on wear of antagonistic and adjacent teeth. *J Adv Prosthodont.* 2020;12(1):1-8. doi:10.4047/jap.2020.12.1.1
46. Makhija SK, Lawson NC, Gilbert GH, et al. Dentist material selection for single-unit crowns: Findings from the National Dental Practice-Based Research Network. *J Dent.* 2016;55:40-47. doi:10.1016/j.jdent.2016.09.010
47. Flanagan D, Ilies H, O'Brien B, McManus A, Larrow B. Jaw bite force measurement device. *J Oral Implantol.* 2012;38(4):361-364. doi:10.1563/AAID-JOI-D-10-00101
48. Javidi H, Tickle M, Aggarwal VR. Repair vs replacement of failed restorations in general dental practice: Factors influencing treatment choices and outcomes. *Br Dent J.* 2015;218(1):E2. doi:10.1038/sj.bdj.2014.1165
49. Raphael KG, Janal MN, Sirois DA, Dubrovsky B. Masticatory muscle sleep background electromyographic activity is elevated in myofascial temporomandibular disorder patients. *J Oral Rehabil.* 2013;40(12):883-891. doi:10.1111/joor.12112
50. Zhou Y, Gao J, Luo L, Wang Y. Does Bruxism Contribute to Dental Implant Failure? A Systematic Review and Meta-Analysis. *Clin Implant Dent Relat Res.* 2016;18(2):410-420. doi:10.1111/cid.12300
51. Baldos J, Szymczak A. DENTIST PRICES and DENTAL COSTS: Finding Affordable Dental Care. DENTIST PRICES and DENTAL COSTS: Finding Affordable Dental Care. Published 2020. Accessed March 16, 2021. <https://www.smilecaredental.ca/blog/dentist-prices-dental-costs-finding-affordable-dental-care>
52. Belle Rive Dental Clinic. Dental Implants Cost. How Much Does a Dental Implant Cost? Published 2020. Accessed March 16, 2021. <https://www.bellerivedental.ca/>
53. Serra CM, Manns AE. Bite force measurements with hard and soft bite surfaces. *J Oral Rehabil.* 2013;40(8):563-568. doi:10.1111/joor.12068
54. Paphangkorakit J, Chaichit R, Khummool K, Noppaisit S, Saengsai A. Effect of biting speed and jaw separation on force used to incise food. *J Oral Rehabil.* 2020;47(6):731-735. doi:10.1111/joor.12961
55. Koc D, Dogan A, Bek B, Yucel M. Effects of increasing the jaw opening on the maximum bite force and electromyographic activities of jaw muscles. *J Dent Sci.* 2012;7(1):14-19. doi:10.1016/j.jds.2012.01.002

